




HEINE HKL® Universal Fiber Optic Examination Light

The ideal light for any doctor's practice. Provides light similar to daylight. Can be powered by any HEINE F.O. projector.


- ∴ **Adjustable Diaphragm Aperture.** Wide range of spot sizes. Bright, even light. No temperature increase.
- ∴ **Colour temperature 3800K.** "Daylight" illumination.
- ∴ **Long, flexible arm "stays put" where required.** Very flexible to position with long-reach arm.
- ∴ **Slimline head.** Unobstructed view.
- ∴ **Universal-mounting.** Suits the practitioner's needs.
- ∴ **Can be retro-mounted to existing projectors.** Makes an examination light out of any projector.
- ∴ **30,000 Lux light output at 20 cm with projector HK7000.**
- ∴ **Mounting options:**
 - on wheeled stand [01]
 - with clamp for table-top mounting (maximum 40mm thick) [02]
 - with clamp for vertical mounting to a tube up to 40mm dia. [03]
 - with clamp for mounting to square tube up to 40mm dia. [04]
 - wall-mounting for HK7000 projector [05]

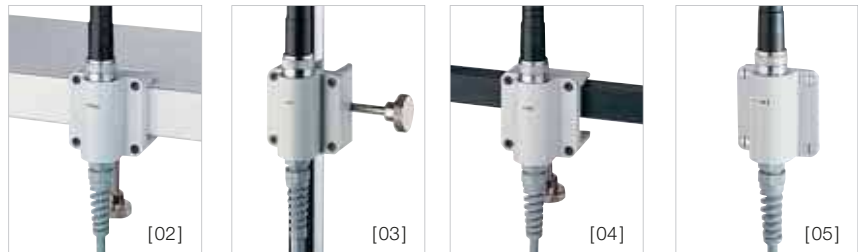

HKL Light	F.O.
on wheeled stand [01]	J-003.27.004
with clamp [02] [03] [04]	J-003.27.001
with wall-bracket [05]	J-003.27.002
Clamp black, only	J-000.27.010
Clamp white, only	J-000.27.012
Base plate, only	J-000.27.011
Wheeled stand without base-plate	Y-096.50.001
Wheeled stand with base-plate	Y-096.50.002

Please order the projector separately. See pages 130–131.

Light spot sizes

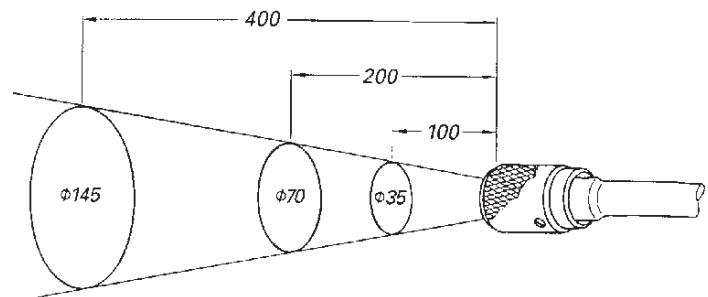

230V is the standard mains voltage. If a different voltage is required, e. g. 120V, please specify when ordering.